TRENDY 971



ready-to-run







medical beds worldwide

CORPORATE PROFILE

- COMPANY IMO-Indústrias Metalúrgicas, S.A.
- > HEAD-QUARTERS Longra, Felgueiras, Portugal
- > VAT No. 500.136.939
- > FOUNDATION YEAR 1946
- > EMPLOYEES 112
- SHARE CAPITAL 2,5 M€
- TURN-OVER (2015) 5,7 M€
- () INSTALLED PRODUCTS 46 países
- NDUSTRIAL FACILITIES 8.400 m2 área coberta
- > PRODUCTION CAPABILITIES Research and development

Prototyping and tools
Molds Engineering
Laboratory quality control
Electronics

Metalworking Plastics

Woodwork

Upholstery
2 assembly line

2 assembly lines

Logistics and distribution

FEATURES



ABS cover base with storage room and O2 bottles compartments.

15

16



⊘ CONTROL **⊘**



COLUMNS ELEVATION SYSTEM Lifting system by 2 hydraulic columns with bilateral pedals control.



TRENDELENBURG
Trendelenburg and AntiTrendelenburg positions by 2
hydraulic columns with bilateral
pedals control.



BACK-REST ARTICULATION Articulated back-rest section by 2 gas springs with 300N capacity with activation rim.





4 SECTIONS ARTICULATIONS Articulated back-rest section by 2 gas springs with 300N capacity with activation rim.

Manually leg-rest articulation by a multipoint mechanical rack.

⊗ ERGONOMY **⊗**



RADIO TRANSLUCENT PLATFORM HPL Compact mattress platform with radio translucent capabilities.

 $2\,or\,4\,sections\,platform\,option.$



X-RAY CASSETTE HOLDER

As an extra option stretcher platform may include an x-ray cassette holder featured to sliding along the mattress platform for bilateral knobs.

It is extractable by the headside to cassette exchange and it has a graduated ruller to adjust on several cassette sizes.



OS DLUX TYPE OPTION

DLUX OPTION

Footboard equipped as an extra option with:

- Foldable monitor tray made in polyethylene.
- HPL Compact footboard customizable color.
- Clinic files tray made in epoxy caoted steel sheet.

⊗ SAFETY **⊗**



CLEANING ENDURANT

The metal frame finished by electrostatic epoxy painting with previous anti-corrosion treatment allows a product suitable for cleaning and hospital disinfecting with alcohol based detergents and moisturized cloth with tepid water.



THE SIDE-RAILS

The stretcher includes a side-rail pair made in epoxy painted steel tube and ABS plastic shoulder extension. Collapsible, with safety lock and grip to handling.



LATERAL FENDERS

TRENDY stretcher sides are equipped by lateral metal fenders soft PVC coated to protect against side impacts.




⊘ MOBILITY **⊘**



Ø200MM CASTERS

High quality Ø200mm casters made by Tente German company. Swivel effect, provided with ball bearings and ABS shielded.

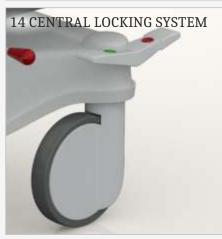

CENTRAL LOCKING

Is provided by 4 brake pedals, one in each bed's corner, with central effect and the following functions:

- Free course;
- Total locking;
- -Steering lock.



5.TH WHEEL OPTION

The fifth wheel system significantly improves the bed's maneuverability in narrow spaces.

The control is released by bilateral pedal.



ABS COVER BASE

Moulded ABS cover to protect base mechanisms and to facilitate the cleaning tasks.

It has storage compartments to patient's personal belongs and 2 rooms for oxygen tanks.



DIMENSIONS (±5%)

Mattress platform	1900 x 650 mm			
External dimensions	2082 x 832 mm			
Platform sections lenght				
Back-rest	678 mm			
Pelvis	390 mm			
Leg-rest	300 mm			
Foot-rest	532 mm			
Minimum/Maximum heigh	565 / 885 mm			
Elevation stroke	320 mm			

ARTICULATIONS (±5%)


Back-rest maximum angle	75º			
Leg-rest maximum angle	40º			
Foot-rest maximum angle	30º			
Cardiac Chair max. angle	90º			
Maximum tilting				
Trendelenburg	15º			
Anti-Trendelenburg	15º			

CAPABILITIES (±5%)

Patient maximum load	350 Kg
Overall bed weight	124 Kg
Elevation Hidraulic Power	4500 N x 2
Back-rest gas lift power	300 N x 2

CONFORMITÉ ET ENREGISTREMENT

Enregistrement officiel comme des Dispositifs Médicaux Classe I et marcation CE. Selon la directive 2007/47/EE de la Union Européenne, annexe IX. Est un produit contrôlé par le système de gestion de la qualité EN ISO 9001:2008.

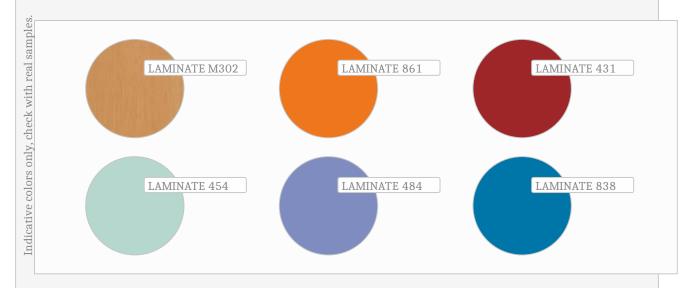

IMO-INDUSTRIAS METALURGICAS, S.A.

www.imo.com.pt

IMO PORTO

Rua Adolfo Casais Monteiro, 117 4050-520 Porto Tel.: 222 009 767

Fax: 222 000 647 porto@imo.com.pt

HEADQUARTERS & FACTORY

Rua Dr. Aurélio Teixeira de Sousa, 575 4650-312 Rande PORTUGAL

T. +351 255 340 220 F. +351 255 340 239

sede@imo.com.pt AUTHORISED DISTRIBUTOR:

IMO LISBOA

Av. Defensores de Chaves, 87-A 1000-115 Lisboa Tel.: 217 963 478 Fax.: 217 975 822 lisboa@imo.com.pt